



Infusion™

Revealing healthy, youthful skin

Incorporating Infusion™ into your clinic opens up a world of opportunity as the treatments can be used in conjunction with other treatment procedures. The perfect complement to any aesthetic procedure, Viora has developed numerous protocols that give you the capability to elevate your offerings, by combining Infusion with Viora's other non-invasive devices and technology.

This advantage allows you to promote the most effective, flexible solutions available. The optimal addition for any aesthetic practice, Infusion™ offers a comfortable and convenient solution for all your patient's skincare needs.

MK-017-ENG

About Viora

Viora, a leading provider of medical aesthetic solutions, is passionate about helping customers deliver the most effective results to their patients. Viora improves upon treatment experiences by enhancing existing technology to give customers the best non-invasive solutions available. With a presence in over 65 countries, Viora's non-invasive products address a wide range of symptoms, treating today's most commonly requested applications for anti-aging, body & facial contouring, cellulite reduction, acne clearance, stretch mark reduction, hair removal, skin lightening and more.

www.vioramed.com | info@vioramed.com |  [vioramed](https://www.facebook.com/vioramed)

**Manufacturer's Address
Viora Inc.**
392 Summit Ave.
Jersey City, NJ 07306 USA
Tel: +1-201-332-4100
Fax: +1-201-332-4200

**Manufacturer's Authorized Representative
Viora UK Ltd.**
450 Bath Road, Longford
Heathrow UB7 0EB, UK
Tel: +44 (0) 208 757 8750
Fax: +44 (0) 208 757 8767



Infusion™

The new wave
of needle-free
mesotherapy

Introducing Infusion™ for Electro-Mesotherapy Trans-dermal Delivery

Getting an edge up on your competitors means that you need to deliver the widest range of treatments and achieve the best results, all for the most attractive pricing. Your offerings need to meet your customers' desires for effectiveness, speed, no pain and no downtime.

Including Infusion™, Viora's electro-mesotherapy device in your clinic provides you with a non-invasive, needle-free method for superior trans-dermal delivery of topicals. Your clients will enjoy the benefits of safe, immediately visible results and revitalized and refreshed skin.

Infusion™ answers your requirements and ensures the optimal delivery of Viora's specially developed topical formulas. Maximum cellular penetration allows you to realize multiple facial and body skin care applications for the best results, every session.

You'll be able to treat the following applications with Infusion™:

- **Anti aging and wrinkle reduction**
- **Collagen regeneration**
- **Skin rehydration**
- **Skin firming & toning**
- **Body sculpting & cellulite reduction**
- **Stretch mark reduction**
- **Skin lightening, anti-aging & resurfacing**
- **Blemish and oil control**



Infusion™ Ampoules - Sterile, High-Quality, Hypoallergenic and Non-Irritating

All natural, Viora solutions are created from the purest, high quality ingredients to maximize cell penetration and ensure safety and effectiveness. Scientifically engineered to achieve the desired effect, each formula is a blend of special serums, proteins and healing agents and is precisely formulated for the specific application. The use of these unique Infusion™ ampoules in conjunction with the Infusion™ system is highly recommended to attain optimal results.

Viora's ampoules include:

Infusion Face Collagen Booster for sagging skin to increase skin's elasticity

Infusion Face Anti-Wrinkling & Lifting specially formulated for mature skin

Infusion Face Intensive Hydration for sensitive and dry skin

Infusion Bright to lighten and improve skin tone

Infusion Smooth for targeted acne treatment

Infusion Body Lipoelim to reduce cellulite

Infusion Body Skin Firm to tighten and replenish dry, sagging skin

Infusion Body Anti-Stretch Mark to treat damaged skin and restore elasticity



"On the advice of my doctor, Dr. Patricia Rittes, I had Infusion treatments done for my face - my skin had been oily, with clogged pores leading to breakouts. Infusion was a simple, painless treatment and now my skin is much clearer, less oily and complete rehydrated. After my sessions with Infusion I can see a big difference - it feels like I have completely new skin!"


- R.P., patient

Technology Infused with Innovation

Infusion™ ionic solutions are driven into the skin by way of electric pulse waves, powered by proprietary Ionwave™ technology. Engineered to work with the skin's biological processes, the Ionwave™ carefully opens micro water channels in the cell membrane to allow the safe passage of the topicals into the cells, bypassing the stratum corneum layer, without damaging the skin.

Ionwave™ Effects:

1. When the Ionwave™ properties come in contact with the cell they begin the structural rearrangement of the lipid bilayer to form micro pores in the cell membrane.
2. As the gatekeepers to the cell, micro channels enable or restrict the passage of natural nutritional compounds and extracts. Micro channels are short-lived and reseal soon after compounds are delivered

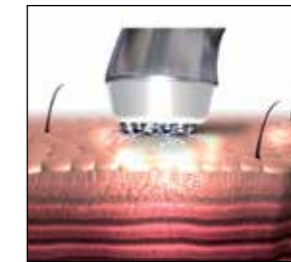

Viora Support

Viora offers the best marketing, clinical and technical support available. These ongoing tools and techniques allow you to be secure in the knowledge that a stable, committed company is partnering with you, for your success.

Technical Specifications

Energy Source	Bipolar balanced output (proprietary wave form)
Output Impedance	10 kilo Ohm
Outputs	2 independent outputs: 1 Face treatment hand piece 1 Body Treatment hand piece
Power Output Adjustment	Linear, continues adjustment from 0 to maximum
Electrical	220 V, 50 c/s or 110 V, 60 c/s-50 V/A
Body Probe Diameter	60mm
Facial Probe Diameter	28mm
Dimensions	35cm (length) x 27cm (width) x 24cm (height)
Weight	5.0 kg /11 lbs

