


Universal IPL

With Optimal Pulse Technology - OPT™

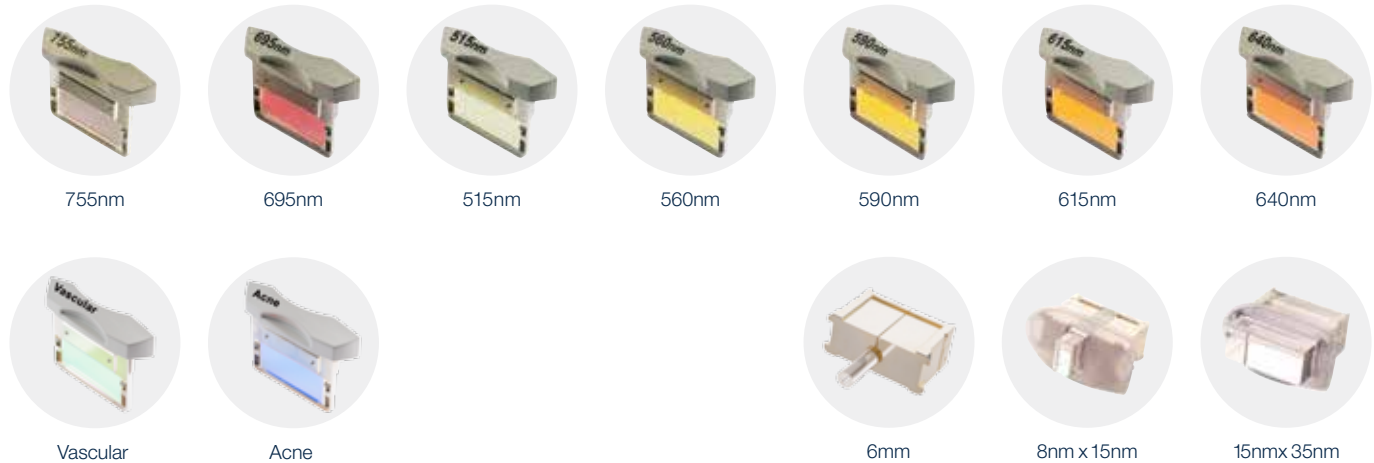


One Applicator for All IPL Indications

9 Filters and 3 Lightguides changeable in seconds

The Lumenis OPT™ advantages

- Control pulse shape by equalizing energy distribution over the entire pulse or sequential pulses
- Gentler, more comfortable, patient-friendly procedures, with lower effective fluences
- High peak power, shorter pulses – optimal for IPL skin treatments using photorejuvenation and treatment of benign pigmented lesions
- Advanced OPT will also allow determining the specific fluence per sub-pulse when using multiple-sequential pulsing, for fine-tuned treatment settings.

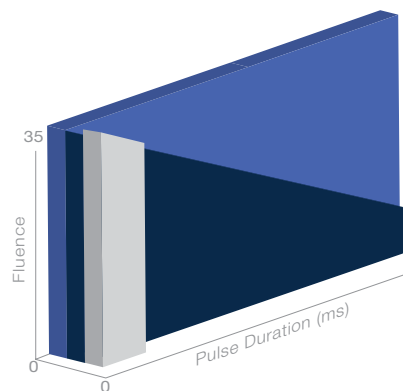


M22 IPL with OPT
High-powered IPL without OPT
Limited capabilities of lower-powered systems


Treat more than 20 skin conditions

Including IPL skin treatments, vascular lesions, pigmented lesions, active acne and hair removal

Typical Clinical Results

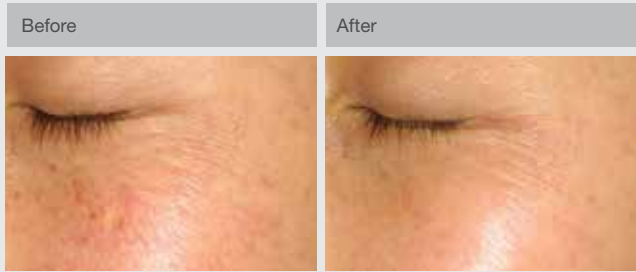

Courtesy of Hiroki Kanda, MD

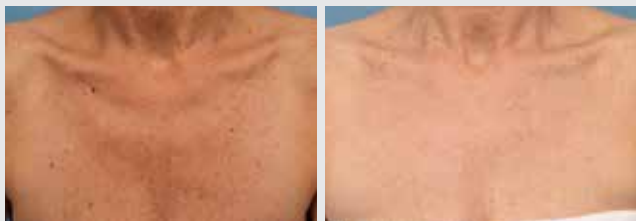

Courtesy of Mellanie Palm, MD



Courtesy of Gilly Munavalli, MD



Courtesy of J. Matthew Knight, MD.



Courtesy of Mariela Nazar, MD.



Courtesy of Mary Lupo, MD.

Warnings and risks: M22™ Universal IPL is contraindicated for patients with active infection, viral, fungal or bacterial diseases. The use of M22 IPL could cause redness, swelling, scarring, change of pigmentation. See the system user manual for a complete list of contraindications and risks.

Specifications

Universal IPL	
Spectrum	400-1200 nm
Expert Filters	515, 560, 590, 615, 640, 695, 755 nm, Acne (400-600 & 800-1200 nm), Vascular (530-650 & 900-1200 nm)
Fluence	10-56 J/cm ²

Pulse Duration	4-20 ms
Pulse Delay	5-150 ms
Pulse Characteristic	Multiple Sequential Pulsing
Repetition Rate	Up to 1 Hz
Spot Size	35 x 15 mm ² ; 15 x 8 mm ² , 6mm
Cooling	Continuous contact cooling

Lumenis Ltd.



Yokneam Industrial Park
Hakidma Street 6, P.O.B. # 240
Yokneam 2069204, ISRAEL
Tel: +972 4 959 9000



Lumenis (Germany) GmbH



Heinrich-Hertz-Str 3 D-63303
Dreieich-Dreieichenhain
GERMANY
Tel: +4 9 (0) 6103 8335 0

Rx ONLY

AMERICAS

San Jose, CA, USA
T +1 408 764 3000
+1 877 586 3647
F +1 408 764 3999

JAPAN

Tokyo, Japan
T +81 3 6743 8300
F +81 3 6743 8301

EMEA

Dreieich Dreieichenhain,
Germany
T +49 6103 8335 0
F +49 6103 8335 300
Roma (RM), Italy
T +39 06 90 75 230
F +39 06 90 75 269
Hertfordshire, UK
T +44 20 8736 4110
F +44 20 8736 4119

ASIA / PACIFIC

Beijing, China
T + 86 10 5737 6677
Gurgaon, India
T +91 124 422 07 95
Kowloon, Hong Kong
T + 852 217 428 00
F + 852 272 251 51

For more information: www.lumenis.com

© 2017 All Rights Reserved. The Lumenis Group of Companies. PB-1003628 Rev D. Lumenis, its logo, M22, OPT, Multiple-Sequential Pulsing, ExpertFilters, SapphireCool lightguides, are trademarks or registered trademarks of the Lumenis Group of Companies. Specifications are subject to change without notice Lumenis is ISO 13485:2012 certified. M22 is CE approved. | M22 is FDA cleared.



LUMENIS.COM/AESTHETIC