

ResurFX™


For True Fractional Non-Ablative
Skin Resurfacing



Proprietary CoolScan™
scanner for customizable
treatments and optimal
results

The only true fractional non-ablative technology


- Extremely versatile state-of-the art CoolScan scanner, for customizable treatments
- Patent pending algorithm for controlled non-sequential pulsing ensures precise, homogenous and safe energy delivery for reproducible results
- With up to 20% coverage in one pass
- Cost effective upgrade with NO disposables

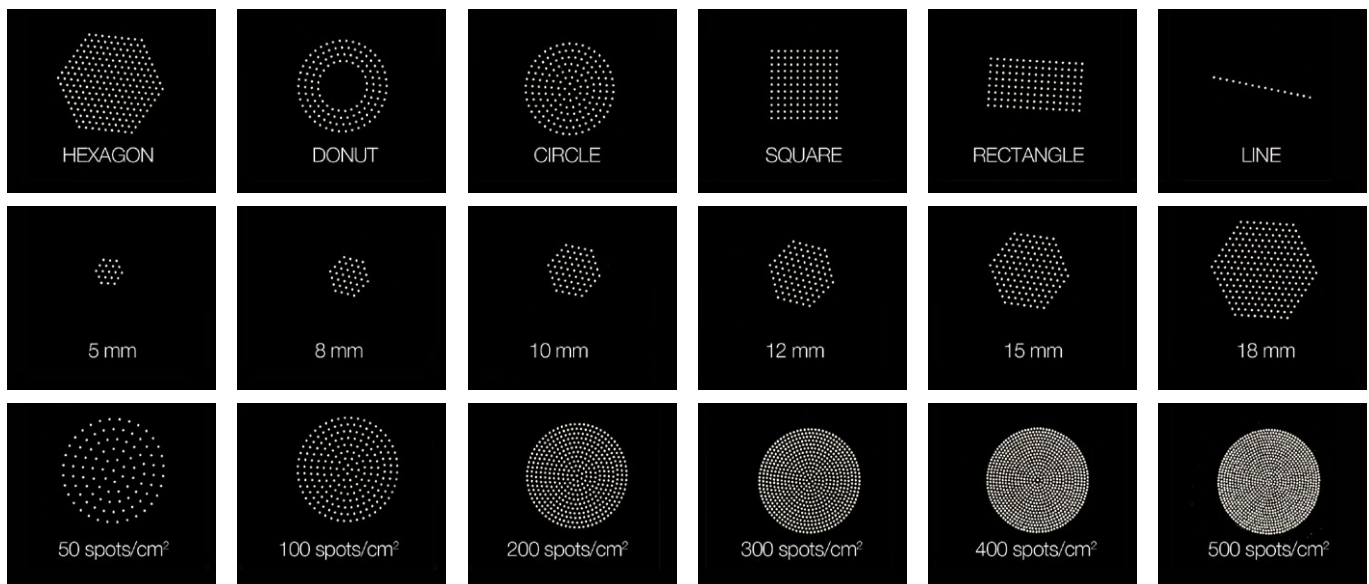


Warnings and risks: ResurFX is contraindicated for patients with any concurrent cancer or history of skin cancer, active infection or bacterial diseases. The use of ResurFX could cause redness, swelling, scarring, and pigmentation change. See the system user manual for a complete list of contraindications and risks.

Fractional non-ablative resurfacing treatments in a single pass

Typical Clinical Results

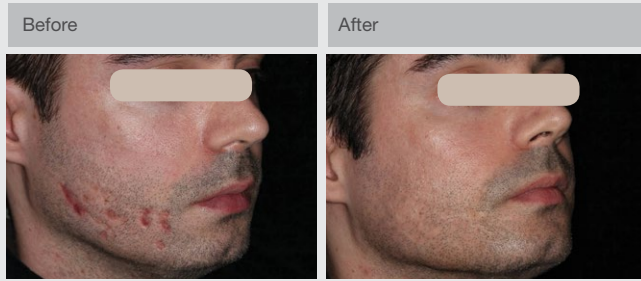

Courtesy of Matteo Tretti Clementoni, MD.

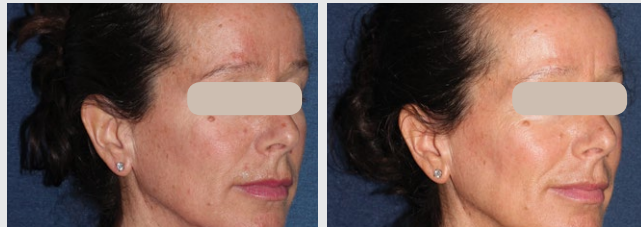

Courtesy of Mitchel P. Goldman, MD



Courtesy of Matteo Tretti Clementoni, MD.



Courtesy of Matteo Tretti Clementoni, MD.



Courtesy of Arielle N.B. Kauvar, MD



Courtesy of Arielle N.B. Kauvar, MD

Specifications

ResurFX™	
Wavelength	1565 nm
Energy	Up to 70 mJ / microbeam
Scan size	Up to 18 mm diameter
Tips	Sapphire tip; Precision tip

Scan shapes	Square, Circle, Rectangle, Line, Hexagon, Donut
Beam density	Up to 500 microbeams / cm ²
Repetition Rate	Up to 2 Hz
Cooling	Continuous contact cooling



Lumenis Ltd.

Yokneam Industrial Park
Hakidma Street 6, P.O.B. # 240
Yokneam 2069204, ISRAEL
Tel: +972 4 959 9000



Lumenis (Germany) GmbH

Heinrich-Hertz-Str 3 D-63303
Dreieich-Dreieichenhain
GERMANY
Tel: +4 9 (0) 6103 8335 0

AMERICAS

San Jose, CA, USA
T +1 408 764 3000
+1 877 586 3647
F +1 408 764 3999

JAPAN

Tokyo, Japan
T +81 3 6743 8300
F +81 3 6743 8301

EMEA

Dreieich Dreieichenhain,
Germany
T +49 6103 8335 0
F +49 6103 8335 300

Roma (RM), Italy
T +39 06 90 75 230
F +39 06 90 75 269

Hertfordshire, UK
T +44 20 8736 4110
F +44 20 8736 4119

ASIA / PACIFIC

Beijing, China
T + 86 10 5737 6677

Gurgaon, India
T + 91 124 422 07 95

Kowloon, Hong Kong
T + 852 217 428 00
F + 852 272 251 51



Rx ONLY



For more information: www.lumenis.com

© 2017 All Rights Reserved. The Lumenis Group of Companies. PB-1000139US Rev D. Lumenis, its logo, M22, ResurFX and CoolScan are trademarks or registered trademarks of the Lumenis Group of Companies. Specifications are subject to change without notice Lumenis is ISO 13485:2012 certified. | M22 is CE approved. | M22 is FDA cleared.



LUMENIS.COM/AESTHETIC